



An ISO 9001:2015
CERTIFIED COMPANY

WORLD SKILL COUNCIL



AUTHORISED
CORPORATE SKILL INSTITUTE

Prife International Pvt Ltd - Brief Introduction

Prife International Pvt Ltd, incorporated in September 2021, stands as the pioneer in bringing innovative Terahertz technology-based devices to India. As the headquarters based in Chennai, the company has established itself as a trusted leader in the wellness and health industry, offering cutting-edge Terahertz products across the nation.

Our diverse range of products includes the Classic, Premium, Power Pro, Magic Patch, and Foot Massager, all powered by Terahertz frequencies, which are revolutionizing the approach to human pain relief and wellness. Prife International is proud to hold a registered trademark, ensuring authenticity and quality in everything we offer, while cautioning customers against counterfeit products in the market.

Beyond products, Prife International is deeply committed to building awareness and educating communities about the significance of Terahertz frequencies. In just a few years, we have built a robust presence and celebrated remarkable achievements. Our success is rooted in our dedication to innovation, customer satisfaction, and fostering a community of wellness advocates.

Prife International Pvt Ltd remains steadfast in its vision to make Terahertz technology a household name, driving wellness and positivity for generations to come. To learn more, visit us at www.iteracare.in.

About the Founders

- **Mr. Bala Mohan Shanker**, the esteemed founder of Prife International Pvt Ltd and iTeracare Pvt Ltd, was a towering figure in the network marketing industry. With over 27 years of full-time experience, he earned a reputation for his integrity, leadership, and unparalleled contributions to the field.
- With a visionary mindset, Mr. Shanker introduced unique products like iTeraCare, aiming to make it a household name globally, particularly in India.
- Tragically, Mr. Bala Mohan Shanker is no longer with us, leaving behind a legacy of dedication, innovation, and inspiration. His vision and values continue to guide Prife International Pvt Ltd, inspiring the organization to carry forward his mission of revolutionizing wellness through Terahertz technology.
- **Mrs. Vedavathi Gorla**, Chairman & Managing Director of Prife International Pvt Ltd, is a dynamic leader with a wealth of experience in business administration and management. With an MBA in Marketing as her academic foundation, she has played a pivotal role in shaping the company's vision and driving its growth as a trusted leader in the wellness and Terahertz technology industry.
- As a seasoned professional, Madam Vedavathi Gorla brings a strategic perspective and a passion for excellence to her role, ensuring that Prife International Pvt Ltd continues to deliver groundbreaking solutions while maintaining its position as a registered trademark and a symbol of trust. Her vision for Prife International extends beyond business growth, focusing on creating lasting value for communities and building a legacy of wellness that resonates across generations.



CERTIFICATE OF REGISTRATION

This is to certify that the Management System of

Prife International Private Limited

New No:9, Old No.23, Plot No.4363, X Block, 5th Street, Anna Nagar West, Chennai - 600040, Tamil Nadu, India.

has been assessed and registered by Lexus Certifications Intl as conforming to the requirements of

ISO 9001:2015

Quality Management System

The Quality Management System is applicable to:

Wellness and Pain Management Technology

Certificate No : LXCI-77-1233

Initial Certificate Date : 18-09-2024
Certificate expiry date : 17-09-2027
1st Surveillance due date : 17-09-2025
2nd Surveillance due date : 17-09-2026

Note: Validity of this certificate is subject to successful completion of each annual surveillance audit on or before due date, in case surveillance audit is not conducted then this certificate shall be suspended/cancelled.

Accredited Board:
Accred Services International Board (Europe) Limited (ASIB)
128 City Road, London, United Kingdom, EC1V 2NX.




Authorised Signatory
Lexus Certifications Intl

This certificate remains valid while the holder maintains the management system in accordance with the standard(s) above, which will be periodically audited by Lexus Certifications Intl.

This certificate remains the property of Lexus Certifications Intl and must be returned on request. In the issuance of this certificate, Lexus Certifications Intl no liability to any party other than to the client, and then only in accordance with the agreed upon certification agreement. Validity of this certificate may be confirmed at www.lexuscertifications.com or email to admin@lexuscertifications.com.

LEXUS

Prife International Pvt. Ltd
has been assessed and registered by

Lexus Certification Intl

confirms the requirements of an

ISO 9001:2015 Company

Quality Management System

applicable to

Wellness and Pain Management Technology



Prife International Pvt. Ltd
has been assessed and certified by
WORLD SKILL COUNCIL
as an authorised
Corporate Skill Institute



WORLD SKILL COUNCIL



WSC, LONDON, UNITED KINGDOM

CENTRAL BHARAT SEVAK SAMAJ

NATIONAL DEVELOPMENT AGENCY, ESTABLISHED IN 1952 BY PLANNING COMMISSION, GOVERNMENT OF INDIA

AUTHORISED CORPORATE SKILL INSTITUTE

This is to certify that

WSC PRIFE INTERNATIONAL PVT LTD
NEW NO.09, OLD NO.23, PLOT NO.4363, X.BLOCK, 5TH STREET, ANNA NAGAR WEST, CHENNAI-
600040, TAMILNADU, INDIA

is the authorised

CORPORATE SKILL INSTITUTE

of the

WORLD SKILL COUNCIL

for the period from 11/Jun/2025 *to* 31/May/2026

Skills affiliated are

FREQUENCY BASED THERAPY, OVER ALL BODY WELLNESS, PAIN MANAGEMENT, ACTIVATION OF
MEREDIAN POINTS, WIEGHT MANAGEMENT, STRESS MANAGEMENT, ANTI AGINF THERAPY,
ANXIETY, FOOT REFLOXOLOGY

Approval Number: WSCCSI-1123

Issued on: 11/Jun/2025



Note: This is a digital certificate. It does not require any signature. To verify authorization details, Scan QR Code.

Started at the instance of Planning Commission, Government of India with the object of enabling individual citizen to contribute, in the form of an organized co-operative effort, to the implementation of the National Development Plan. The constitution and functioning of Bharat Sevak Samaj is approved unanimously by the Indian Parliament.

www.worldskillcouncil.org



Office Location & Contact Details

Corporate Office:

**New No.09, Old No.23,
Plot No.4363, X Block,5th Street,
Anna Nagar West,
Chennai-600040**

**Phone No: 044 – 2626 7777
78458 76787 | 9940036787**

**www.iteracare.in
prifeinternational777@gmail.com**







Regional Office:

**Sangam Apartment
No: 1539/4, J Block, 10th Street,
Anna Nagar West,
Chennai-600040**








**Phone No: 044 – 2626 7777
78458 76787 | 9940036787**

**www.iteracare.in
prifeinternational777@gmail.com**

The Future of Wellness: Terahertz and Quantum Therapy

-  **Terahertz frequencies are found between far-infrared rays and microwave energy. It refers to electromagnetic waves with frequencies in the range 0.1 to 10 THZ.**
-  **Terahertz (Thz) light waves are radio light waves of very short wavelength whose frequency ranges from 0.03 thz to 0.3 thz.**
-  **1 THz to 2 THz: This range is thought to be especially beneficial for healing processes.**
-  **Researchers have been discovering the health benefits of terahertz therapy on the human body.**
-  **Often referred to as "T-rays", terahertz waves exhibit unique properties that make them suitable for a variety of applications, including imaging, sensing, and communication.**
-  **Terahertz waves can benefit 3-5 centimeters on human skin, directly act on the deep tissue, supplement the energy of life cells, regulate the physiological functions of the body and accelerate the generation and decomposition of cells.**

The Future of Wellness: Terahertz and Quantum Therapy

-  **The Terahertz wave brings in a large number of negative ions, solving the problem of excessive accumulation of positive ions, further enhancing the recombination ability.**
-  **It can promote blood circulation, reduce inflammation and relieve pain, and has remarkable curative effect on various inflammation and rheumatism.**
-  **Stimulate the body's cells self-repair and proliferation of the ability to help blood circulation to accelerate the elimination of toxins, from the inside out to solve sub-health problems.**
-  **Terahertz waves act on deep damaged tissues of the human body.**
-  **Improve blood circulation, increases blood flow, eliminates swelling , accelerates tissue healing; reduces muscle tension and improve microcirculation.**
-  **It has good healing effect on bone and joint diseases, muscle and soft tissue sprains, inflammation, neuralgia, traumatic hematoma, etc.**
-  **It can penetrate deep into the subcutaneous tissue and its therapeutic effect can last long after treatment.**

Terahertz Product Pricing and Variants



iTeraCare Classic



iTeraCare Premium



iTeraCare Power Pro



Amazing Body Rejuvenator



5G Magic Patch

Pricing and Variants



iTeraCare Classic

INCREDIBLE CLASSIC 3.0

The terahertz quantum wand is a powerful leap ahead of the terahertz technology.

The Classic Terahertz Device is a revolutionary wellness tool designed to harness Terahertz frequencies, promoting overall health and well-being. It is ideal for personal use and offers the following benefits:

Key Features

Advanced Technology: Utilizes Terahertz frequencies to rejuvenate and energize cells.

Portable Design: Lightweight and easy to use anywhere, anytime.

Affordable Wellness Solution: Cost-effective yet powerful for daily health routines.

User-Friendly: Simple operation with intuitive controls, making it suitable for all ages.

Pricing and Variants

PREMIUM DOMINATOR

The **Terahertz Premium** is a high-performance wellness device designed for those seeking a **more powerful** and refined experience in their health routine. With enhanced terahertz wave emission, it delivers **deeper therapeutic effects**, providing **faster recovery, pain relief, and overall body wellness**. Ideal for users who require a higher intensity for more focused and effective therapy.

The **Terahertz Premium** is designed for users seeking deeper, more efficient therapeutic effects. Whether you need **faster muscle recovery, relief from soreness, or improved sleep**, this advanced device offers customized therapy to help you achieve your **wellness goals** with ease and comfort.



Premium Dominator

Pricing and Variants



iTeraCare Power Pro

ITERACARE POWER PRO

The **Terahertz Pro** is a professional-grade wellness device designed for advanced therapeutic use. Equipped with the latest **terahertz technology**, it provides highly targeted and deep-reaching benefits, making it ideal for **chronic pain management, muscle recovery, and advanced therapeutic treatments**. The Terahertz Pro is a powerful, **advanced device** designed for **professional-level therapy, offering precision, customization, and larger coverage** for comprehensive health and wellness treatments.

Pricing and Variants

AMAZING BODY REJUVENATOR

The Terahertz Foot Massager is specifically designed to target the feet, a key area for reflexology. By using terahertz waves, it helps to enhance circulation, relieve foot pain, and promote overall health and well-being. Ideal for anyone looking to experience relaxation and healing through a therapeutic foot massage.

The Terahertz Foot Massager is the perfect solution for those looking to relieve foot pain, improve circulation, and promote overall health through therapeutic foot care. Whether at home or at work, enjoy the soothing benefits of terahertz wave therapy at your convenience.

Transfer the Micro magnetic vibration energy to the whole body through the foot electrode plate.

The cells produce polarized heat and the positive and negative ions in the body move violently.

Deep Physical Therapy of Human Muscle Energy from Inside out.
Terahertz radiation could potentially penetrate tissues to a certain depth without causing harm.

Terahertz technology could potentially be used to enhance the massage experience by penetrating deeper into the tissues of the feet. This could potentially help in relieving muscle tension, improving circulation, and providing overall relaxation.



Amazing Body Rejuvenator

Pricing and Variants



5G MAGIC PATCH

Targets and alleviates joint pain, muscle soreness, and discomfort.

Compact, lightweight, and adhesive, allowing you to wear it anywhere, anytime, on-the-go relief and reusable up to 200 times.

No chemicals, No side effects—just natural frequencies to support your body’s healing process.

No bulky devices or complicated instructions—just stick it and experience the benefits.

No batteries, chargers, or extra maintenance required.






Works discreetly under clothing for seamless integration into your daily routine.

Understanding the differences

A Breakdown of our Terahertz Products

Features	Classic, Premium & Pro Models	Foot Massager
Target Area	Multi-area body therapy (Back, Neck, Shoulders, etc.,)	Foot-specific therapy etc.,
Design	Handheld or wand-like, for bodywide use	Ergonomically designed for feet
Functionality	General pain relief, relaxation, detoxification, skin rejuvenation	Focuses on reflexology, foot pain relief, and circulation
Customization	Multiple intensity and therapeutic modes for varied body areas	Adjustable settings for foot massage and heat
Portability	Portable, suitable for different areas and on-the-go use	Stationary, Home use only
Therapeutic Focus	Muscle relief, detoxification, skin rejuvenation, circulation	Reflexology, foot pain, and circulation improvement
Best for	Full-body therapy, general wellness, and muscle recovery	Foot care, relaxation, and pain relief

Technical Parameters

Features	Voltage	Frequency	Power	Energy Gear
 Classic	110V/220V	50HTz/60HTz	750 Watts	3 Modes
 Premium	110-130V/ 220-250V	50HTz/60HTz	750 Watts	3 Modes
 Power Pro	220V	50HTz	1200 Watts	3 Modes
 Foot Massager	220V	-	260 Watts	-
 Magic Patch	-	-	-	-

Variants & Prices

iTeraCare Classic



32,450/-

Amazing Body Rejuvenator



3,24,500/-

iTeraCare Power Pro



3,24,500/-

Premium Dominator



97,350/-

5G Magic Patch



13,000/-

Mentioned prices inclusive of GST

Terahertz Technology for Humankind



Cellular-Level Healing and Regeneration Energy

Terahertz waves resonate with the natural frequency of human cells, promoting better cellular communication and repair.

Improved Cellular Function: Encourages detoxification by enhancing the vibration of water molecules, aiding in nutrient absorption and waste elimination.



Stress Reduction and Relaxation Non-Invasive Therapy

Devices using THz frequencies are designed to target specific areas of discomfort, helping to relax tense muscles and reduce chronic stress.

Calming Effect: THz frequencies have been associated with balancing the body's natural bioenergy, fostering a sense of calm and well-being.



Pain Management Without Side Effects Chronic Pain

Terahertz-based devices offer relief for conditions such as arthritis, migraines, and fibromyalgia by reducing inflammation and improving blood flow.

No Medication required: Provides a drug-free alternative for managing pain, reducing dependence on pharmaceuticals and avoiding side effects.

TeraHertz Technology for Humankind



Boosting Immunity and Vitality Immune System Activation

The technology aids in stimulating white blood cells, which enhances the body's natural defence mechanisms.

Energy Revitalization: Encourages better oxygenation and circulation, resulting in improved energy levels and overall vitality.



Skin and Tissue Rejuvenation Anti-Aging Benefits

THz frequencies improve skin elasticity, promote collagen production, and reduce the appearance of wrinkles.

Wound Healing: Speeds up the healing of wounds and injuries by increasing cellular regeneration and reducing infection risks.



Detoxification and Hydration Water Molecule Activation

Enhances the structure and vibrational energy of water, allowing cells to stay hydrated and detoxified effectively.

Toxin Removal: Supports the body's natural detox pathways, leading to better digestion and organ function.

TeraHertz Technology for Humankind



Hormonal Balance and Endocrine Support Regulation of Hormones

THz devices have been observed to help balance hormonal levels, benefiting conditions like thyroid dysfunction or hormonal imbalances during menopause.

Mental Clarity: Reduced hormonal fluctuations lead to better mood stability and focus.



Enhanced Circulation and Oxygenation Improved Blood Flow

Helps in dilating blood vessels and reducing blood viscosity, ensuring better circulation and oxygen delivery to vital organs.

Organ Function Optimization: Better oxygenation results in improved heart, brain and kidney function, promoting overall wellness.



Support for Chronic Conditions Diabetes Management

Encourages improved cellular insulin sensitivity and glucose metabolism.

Cardiovascular Wellness: Reduces inflammation in arteries and promotes better heart health through improved circulation.

Guide for Effective TeraHertz Therapy



Relief for Varicose Veins with Terahertz Therapy


Relief for Varicose Veins with Terahertz Therapy A non-invasive solution for pain, swelling & discomfort.

What Are Varicose Veins?

Varicose veins are enlarged, twisted veins, usually in the legs, caused by weakened vein walls and poor blood flow.

They often cause:

- ⊗ Swelling and heaviness
- ⊗ Leg pain or burning sensation
- ⊗ Visible, bulging veins
- ⊗ Fatigue or cramps, especially after standing

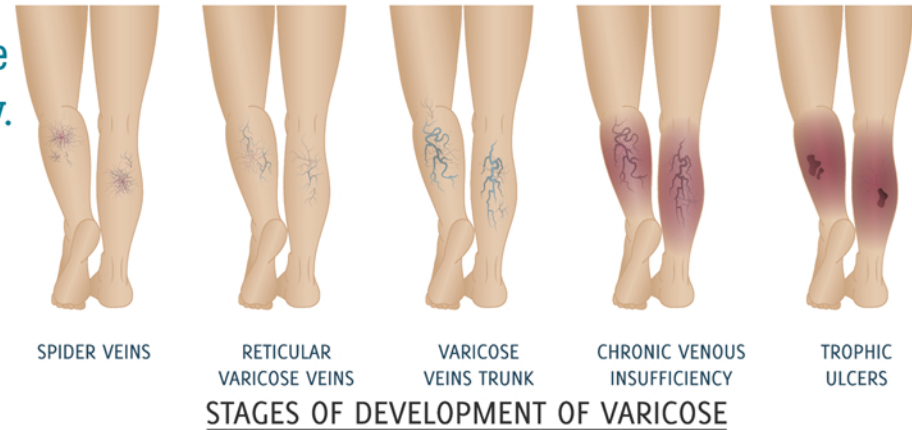

How the Terahertz Device Helps Our Terahertz Frequency Device works on a cellular level using safe, non-ionizing waves and soothing thermal energy to:

1. Improve Blood Circulation Stimulates better venous return and helps reduce pressure build-up in the veins.
2. Reduce Inflammation & Pain Calms the inflamed blood vessels, easing the pain and heaviness in the legs.
3. Relieve Swelling Assists in fluid drainage and lowers swelling around the ankles and calves.
4. Energize Vein walls Enhances cellular activity, promoting stronger and more flexible blood vessels over time.
5. Support Detoxification Removes trapped waste products that contribute to vein damage and discomfort.

Guide for Effective TeraHertz Therapy

Where to Blow the Device for Varicose Veins

Surrounding the Affected Area (Not directly on protruding veins):

- ◆ Gently blow around the varicose veins in circular motion (not on the raised veins directly).
- ◆ Duration: 15–20 mins per leg

Behind the Knees:

- ◆ Stimulates lymphatic and venous drainage
- ◆ Duration: 3–5 mins

Groin Area (Inguinal Nodes):

- ◆ Improves blood outflow from legs
- ◆ Duration: 3–5 mins

Soles of the Feet:

- ◆ Stimulates circulation from the ends of the body
- ◆ Duration: 3 mins per foot

Precautions:

- ⌚ Do not blow directly on bulging veins – always go around them.
- ⌚ Avoid long exposure to high heat on one spot.
- ⌚ Drink charged water before and after use.
- ⌚ Best results with consistent daily use.
- ⌚ Wear compression stockings if advised by a doctor.
- ⌚ Natural. Non-invasive. Powerful.



Guide for Effective TeraHertz Therapy

Where to Blow the Device for Diabetes

Main Areas:

- ◆ Blow on the abdomen, especially around the navel.
- ◆ Focus on the lower back, near the spine (Bladder 23).
- ◆ Target the inner leg, 3 inches above the ankle (Spleen 6).

Additional Tip: Blow on the soles of the feet to improve circulation.

Where to Blow the Device for Kidney Stones

Main Areas:

- ◆ Blow on the lower back (around the kidney area).
- ◆ Focus on the lower abdomen, near the pubic bone (Ren 3).
- ◆ Target the sides of the lower back where the kidneys are located.

Where to Blow the Device for Gall Bladder Stones

Main Areas:

- ◆ Blow on the upper rightside of the abdomen, just below the ribs (liver and gall bladder area).
- ◆ Focus on the back, right side, near the rib cage.



Guide for Effective TeraHertz Therapy

Where to Blow the Device for Parkinsons Disease

Main Areas:

- ◆ Blow on the back of the neck (Gallbladder 20) and the base of the skull.
- ◆ Focus on the hands and feet to improve mobility.
- ◆ Use the device on the spine to relax the nervoussystem.

Where to Blow the Device for Frozen Shoulder

Main Areas:

- ◆ Blow on the shoulder joint directly
- ◆ Focus on the neck and upper back to relax surrounding muscles.
- ◆ Blow along the arm and shoulder blade for improved blood flow.

Where to Blow the Device for Heart Problems

Main Areas:

- ◆ Blow on the center of the chest (Lung 1 area).
- ◆ Focus on the wrist crease near the little finger (Heart 7)
- ◆ Blow gently along the back near the spine, especially the middle-upper back.



Guide for Effective TeraHertz Therapy

Where to Blow the Device for Joint Pain and Arthritis

Main Areas:

- ◆ Blow directly on the painful joint.
- ◆ Focus on nearby areas like below the knee (Stomach 36) and hand (Large Intestine 4).



Where to Blow the Device for Digestive Problems

Main Areas:

- ◆ Blow on the middle of the abdomen, above the navel (Conception Vessel 12)
- ◆ Target the lower abdomen to stimulate digestion

Where to Blow the Device for Stress and Anxiety

Main Areas:

- ◆ Blow on the forehead (between the eyebrows).
- ◆ Focus on the back of the neck and the wrists.



Guide for Effective TeraHertz Therapy



General Relaxation and Calmness

- **Forehead:** Blow gently on the center of the forehead (Yintang, also known as the "Third Eye" point). This can promote calmness and improve focus.
- **Base of the Skull:** Target the area at the back of the neck (Gallbladder 20), which is linked to mental relaxation.



Nervous System Stimulation

- **Spine:** Blow along the entire length of the spine to improve communication within the nervous system. Focus especially on the upper spine (near the neck) and lower back.



Hormonal Imbalances and Menstrual Issues

- **Menstrual Cramps:** Target the lower abdomen (Ren 3) and the lower back (Bladder 23) to reduce pain and improve circulation.
- **Hormonal Balance:** Blow on Spleen 6 (inner leg above ankle) and the abdomen.



Respiratory Problems (Cough, Cold, Asthma)

- **Chest Congestion:** Blow on the chest area, particularly below the collarbone where Lung 1 (LU1) is located.
- **Throat Discomfort:** Direct airflow along the neck and throat area.

Guide for Effective TeraHertz Therapy



General Areas for All Types of Paralysis

- **Spine:** Blow along the entire spine to improve nerve function and circulation. Focus on the affected region (e.g., cervical for upper body paralysis, lumbar for lower body).
- **Base of the Skull (Gallbladder 20):** Target the area at the back of the neck to stimulate the nervous system.
- **Hands and Feet:** Focus on the palms, soles, and the fleshy area between the thumb and index finger (Large Intestine 4).



Hemiplegia (One Side of the Body Paralysis)

- **Affected Side:** Blow directly on the paralyzed limbs, starting from the upper areas (shoulders, hips) and moving downwards.
- **Key Points:** Focus on joints like the shoulder, elbow, wrist, hip, knee, and ankle.



Paraplegia (Lower Body Paralysis)

- **Lower Back:** Blow on the lower back (Bladder 23) and along the lumbar spine.
- **Legs and Feet:** Target the entire length of both legs, especially the soles of the feet and inner legs (Spleen 6).

Guide for Effective TeraHertz Therapy

Quadriplegia (Full Body Paralysis)

- **Spine:** Blow along the entire spine, focusing on areas that control the limbs.
- **Hands and Feet:** Stimulate the palms, soles, and fingers to improve energy flow and circulation.
- **Chest:** Blow on the chest to improve breathing and circulation.(Large Intestine 4).



Facial Paralysis(Bell's Palsy)

- **Face:** Blow gently on the affected side of the face, especially along the jawline, cheeks, and forehead.
- **Neck:** Focus on the sides of the neck to improve circulation to the face.
- **Consistency:** Regular sessions are crucial for effective results.
- **Gentle Treatment:** Use the device gently and avoid direct contact with sensitive skin or open wounds.
- **Patience:** Paralysis recovery takes time, and combining therapy with physical rehabilitation may enhance outcomes.



Usage Guidelines & Responsibilities



General Guidelines for Using Terahertz Devices

- **Preparation:**
 - Ensure the device is clean and operational.
 - The user should be in a comfortable, relaxed position.
 - Hydration: Drink charged water before and after each session to help flush out toxins.
- **Usage Instructions:**
 - Use the device 2-3 times a day for 10–15 minutes per session.
 - Maintain a distance of 2-3 cm from the skin while using.



Disclaimer

- ⚠ iTeracare blower is not a medical device, no medical claims by the manufacturer or on behalf. It is not intended for use in the diagnosis, treatment, mitigation, or prevention of disease.
- ⚠ This product is not recommended for use by pregnant women, children, or women during their menstrual cycle; it may not be suitable for individuals who have undergone cataract surgery or any metal implants in the body or open/bleeding wounds.
- ⚠ Do not use the product/device in the vicinity of stents and pacemakers.
- ⚠ The use of Terahertz devices is intended to support general wellness and is not a substitute for professional medical advice, diagnosis, or treatment. The response to these devices, including the time required for noticeable benefits, may differ significantly among individuals. Variations in outcomes are influenced by factors such as individual health conditions, age, lifestyle, and physiological differences. Prife International Pvt Ltd does not guarantee specific results.

Social Media Handles & Account Details

Please follow us:

-  Prife iTeracare
-  prife_international_pvt_ltd
-  @Prife_Intl777
-  +91 78458 76787
-  @prifeinternationaloffl

BANK DETAILS:

HDFC BANK

BANK A/C NO: 5925 66666 77777

IFSC CODE : HDFC0000017

Warranty Details

A one-year warranty is provided from the date of purchase.

Warranty Coverage: During the one-year warranty period, if the product develops any defects or malfunctions that are not a result of user misuse or damage, the manufacturer or seller will typically repair or replace the product at no cost to the customer.

Proof of Purchase: You will likely need to provide proof of purchase, such as a receipt or invoice, to take advantage of the warranty.

Warranty Exclusions: Warranties usually do not cover damage or issues resulting from misuse, abuse, neglect, accidents, unauthorized repairs or modifications, or normal wear and tear.

Warranty Registration: You may require to register the product to activate the warranty. Be sure to follow any registration instructions if they apply to your product.

Contact Information: Keep the contact information updated in your backend handy in case you need to initiate a warranty claim.

Warranty Period: Make sure you are aware of when the warranty period begins. It typically starts from the date of purchase, but it's always a good idea to confirm this with the authorized distributor.

Maintenance and Care: Proper maintenance and care of the product in accordance with the user manual can help ensure that you don't void the warranty. Neglecting proper care may give the service team a reason to deny a warranty claim.

Warranty Details

Extended Warranty: You may contact our customer care for extended Warranty for an additional cost, which can extend the warranty coverage beyond the standard one year.


Crystal Warranty: "Warranty does not apply to any physical damage or breakage of the crystal of the device."

Warranty will be void under the following conditions:

- ⌚ Exceeded warranty period
- ⌚ Damage due to not following the user's manual
- ⌚ Device damage due to mishandling or repair by unauthorized persons

If your product experiences issues covered by the one-year warranty within the warranty period, you should contact the customer care to initiate a warranty claim. They will provide instructions on how to proceed, whether it's sending the product in for repairs or arranging for a replacement


Autonomic Nervous System

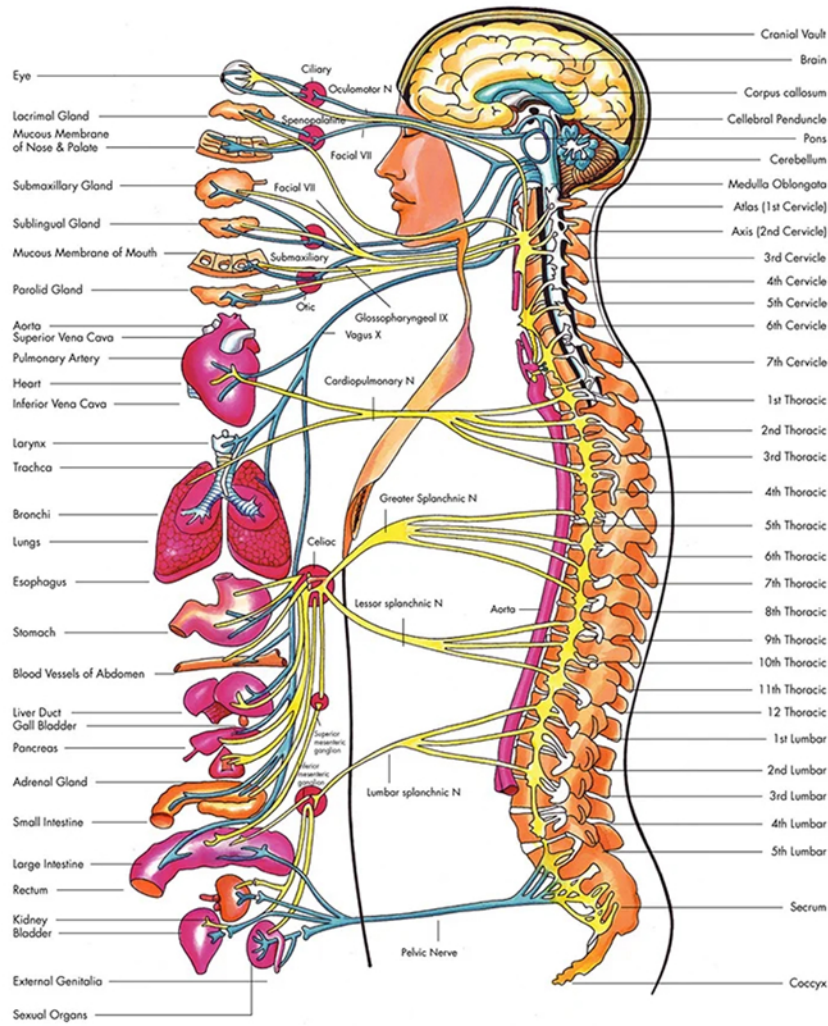


Spinal Nerve Function

VERTEBRAL LEVEL	AREA OR ORGAN	POSSIBLE SYMPTOMS
C1	MASTER CONTROL FOR ENTIRE NERVOUS SYSTEM BRAIN • PITUITARY • FACE EARS • EYES • NOSE • HEAD • PARATHYROID	HEADACHES • MIGRAINES EARSACHES • COLDS/FLU SORE THROAT • SINUSITIS FATIGUE • VERTIGO DIZZINESS • ANXIETY ALLERGIES • MEMORY LOSS VISUAL OR AUDITORY DISTURBANCES
C2	DIAPHRAGM • NECK MUSCLES SCALP • SINUSES • TONGUE TEETH • JAW • THYROID	WEIGHT GAIN • COLD INSOMNIA • FATIGUE COLDS/FLU • SINUSITIS HEADACHES • SORE THROAT VERTIGO • DIZZINESS ALLERGIES • REFLUX NUMBNESS • TINGLING
C3	NECK MUSCLES • THYROID SINUSES • MOUTH VOCAL CORDS • TONSILS	WEIGHT GAIN • COLD INSOMNIA • FATIGUE COLDS/FLU • SINUSITIS HEADACHES • SORE THROAT VERTIGO • DIZZINESS ALLERGIES • REFLUX NUMBNESS • TINGLING
C4	ARMS • ELBOWS • WRISTS HANDS • FINGERS PARATHYROID • ESOPHAGUS ACROSS SHOULDERS	WEIGHT GAIN • COLD INSOMNIA • FATIGUE COLDS/FLU • SINUSITIS HEADACHES • SORE THROAT VERTIGO • DIZZINESS ALLERGIES • REFLUX NUMBNESS • TINGLING
C5	HEART • LUNGS • THYMUS BRONCHIAL TUBES • TRACHEA CHEST • BREAST • PLEURA GALL BLADDER • ARMS ACROSS SHOULDERS UPPER BACK SHOULDER BLADES	HEART CONDITIONS • ASTHMA BRONCHITIS • INDIGESTION DIFFICULTY SWALLOWING NAUSEA • HEADACHES FATIGUE • PLEURISY CHRONIC COUGHING SHORTNESS OF BREATH BLOATING AFTER EATING
C6	STOMACH • PANCREAS SPLEEN • LIVER DUODENUM PERITONEUM • MIDDLE BACK	GASTRITIS • ULCERS INDIGESTION • CRAVE SWEETS HEADACHES • FATIGUE FEELING TOXIC • RASHES CATCHES COLDS EASILY DIABETES • HYPOGLYCEMIA
C7	ADRENALS SMALL INTESTINE KIDNEYS MIDDLE LOWER BACK	STRESS • FATIGUE IRRITABLE BOWELS • HIVES ECZEMA • RASHES HIGH BLOOD PRESSURE GASSY • BLOATING WATER RETENTION
T1	LARGE INTESTINE SMALL INTESTINE ILEOCECAL VALVE APPENDIX • BLADDER OVARIES • UTERUS PROSTATE • TESTICLES REPRODUCTIVE ORGANS RECTUM • URETHRA LOW BACK • PELVIS BUTTOCKS • GROIN THIGHS • KNEES ANKLES • FEET • TOES SCIATIC NERVE	CONSTIPATION • DIARRHEA GASSY • BLOATING IRRITABLE BOWEL BLADDER INFECTIONS PAINFUL URINATION BED WETTING HEMMORHOIDS MENSTRUAL PROBLEMS PREMENSTRUAL SYNDROME IRREGULAR PERIODS MENOPAUSE INFERTILITY LOW SEX DRIVE NUMBNESS OR TINGLING YEAST INFECTIONS SCIATIC PAIN SPRAINED ANKLES
T2		
T3		
T4		
T5		
T6		
T7		
T8		
T9		
T10		
T11		
T12		
L1		
L2		
L3		
L4		
L5		
SACRUM		



Ear Reflexology Chart

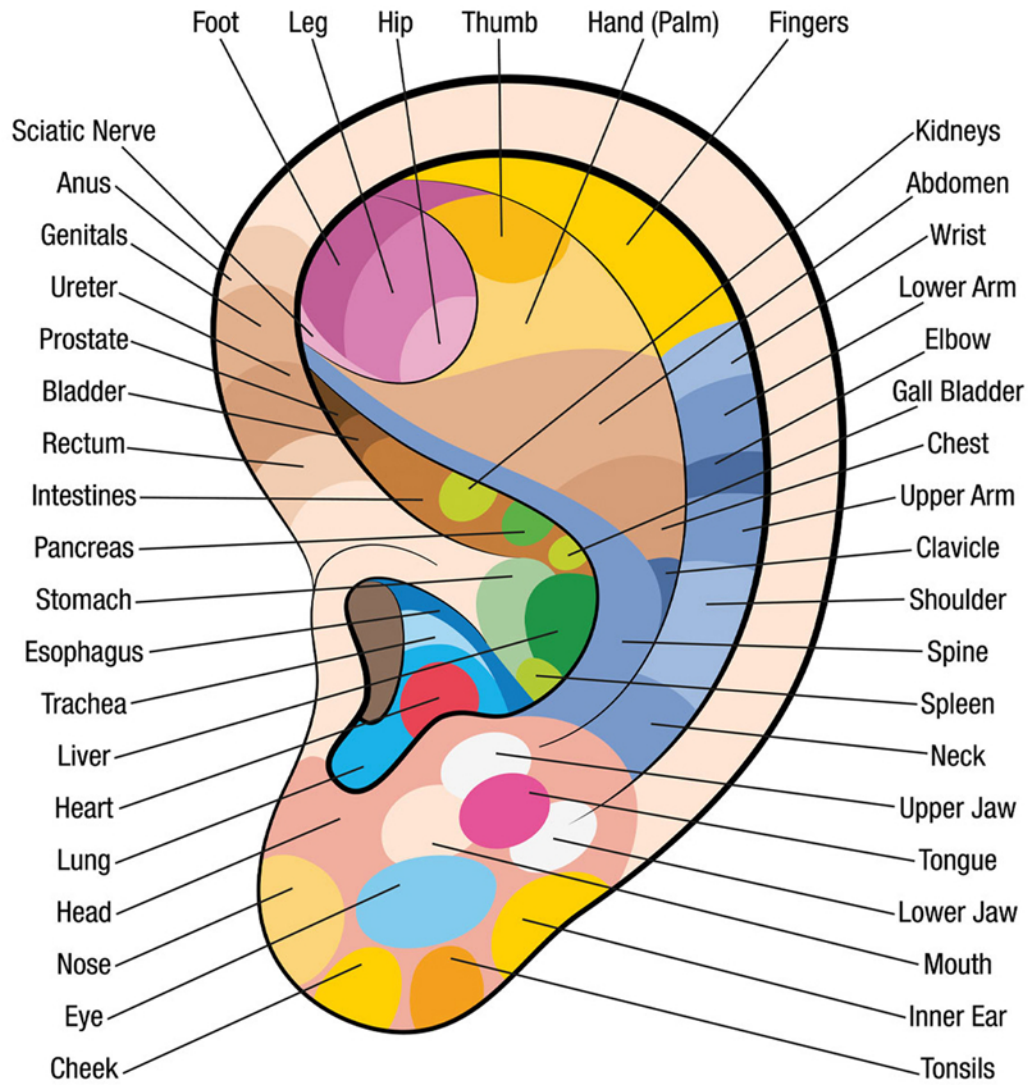





Hand Reflexology Chart

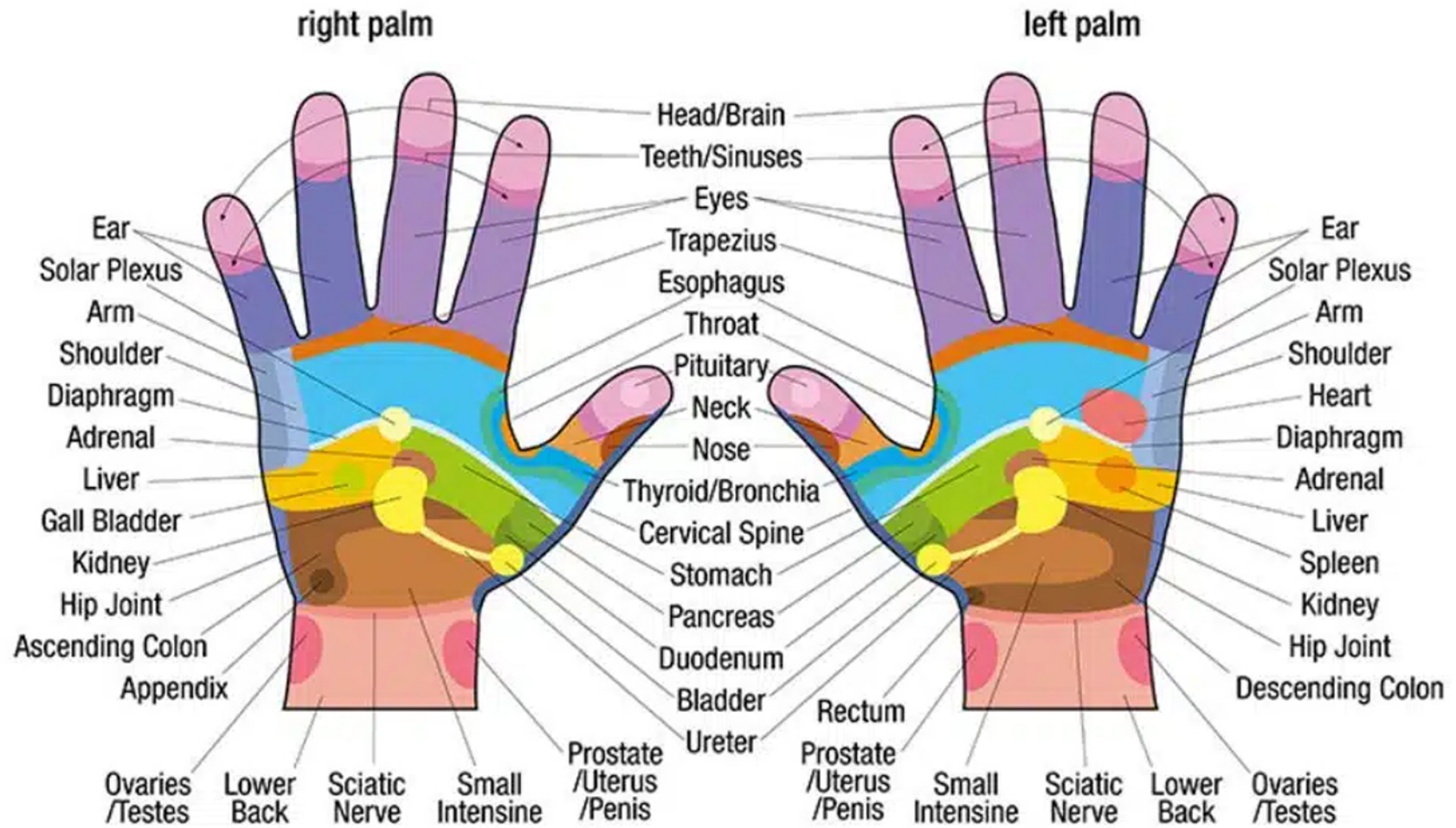

Foot Reflexology Chart

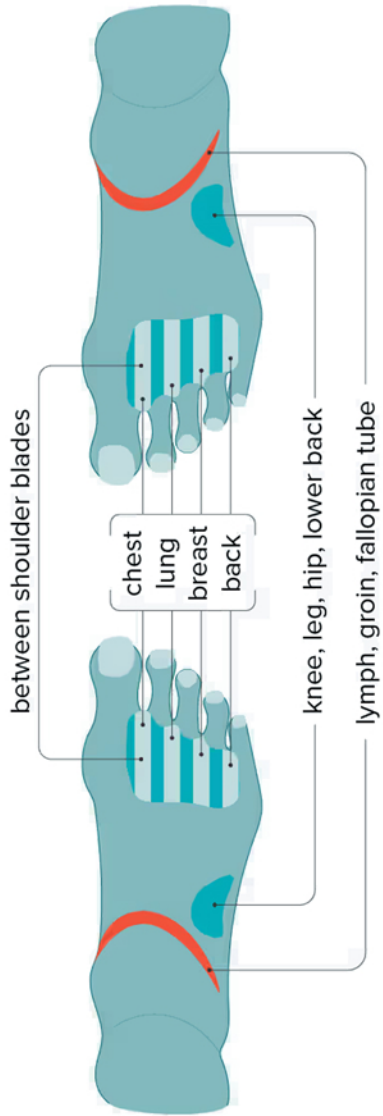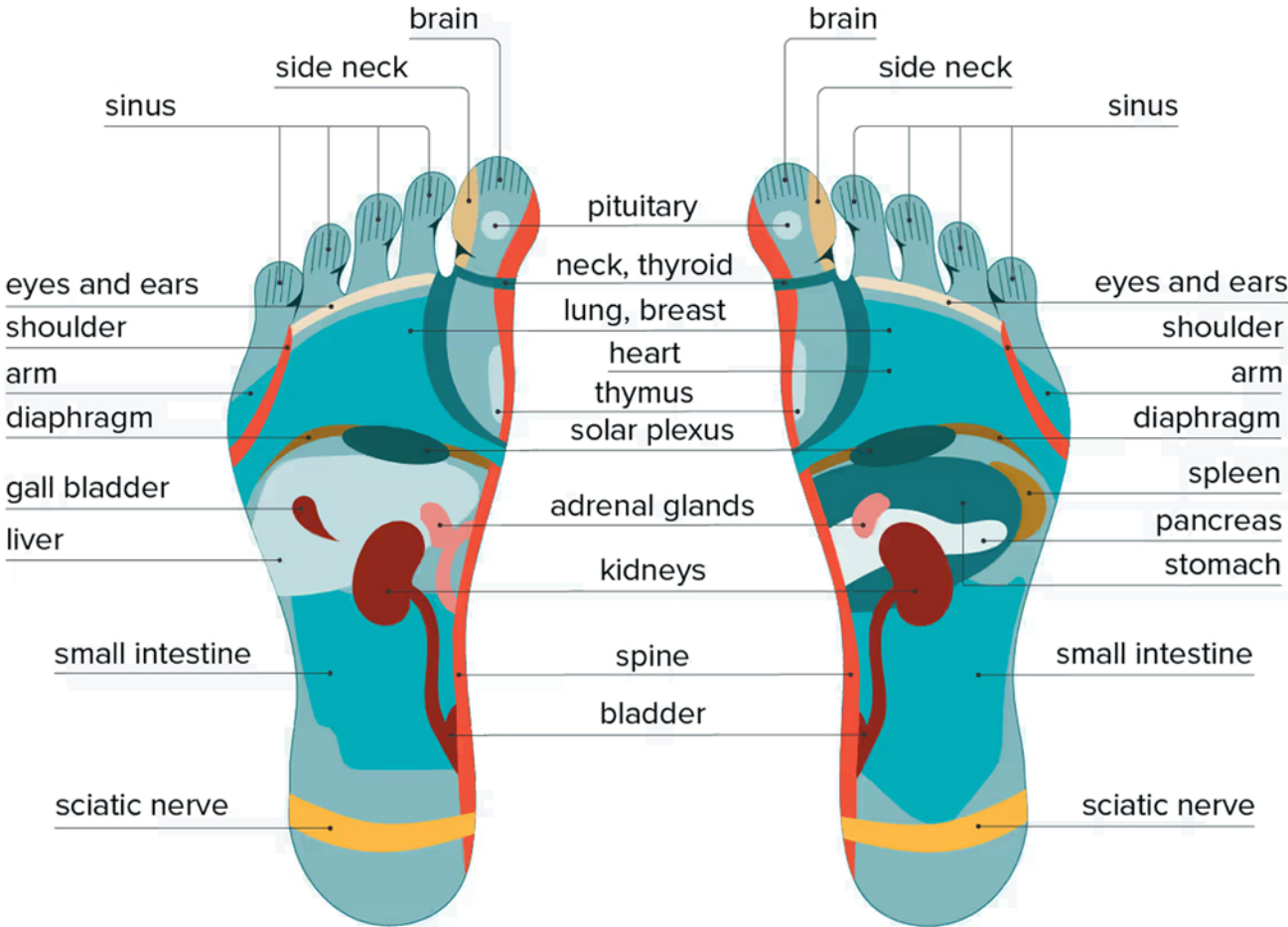


Foot Reflexology Chart

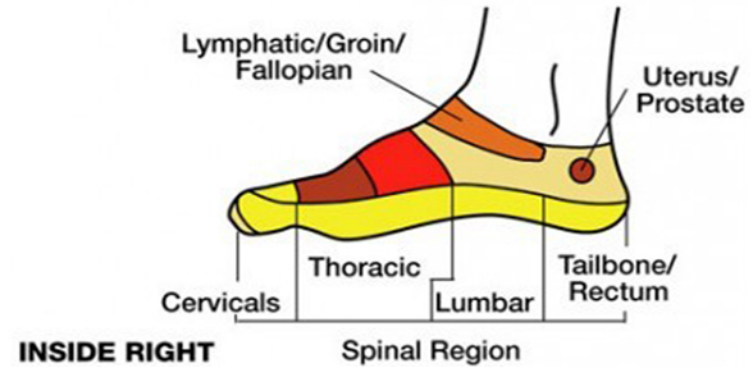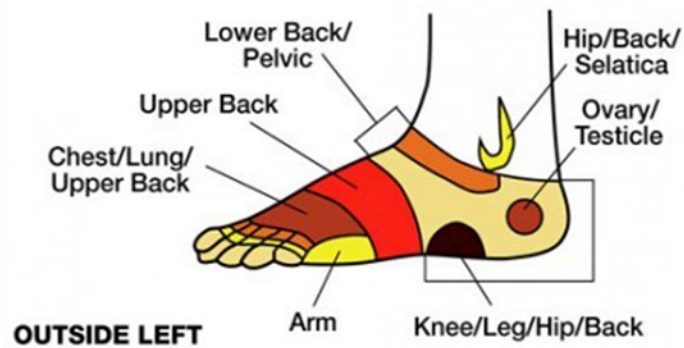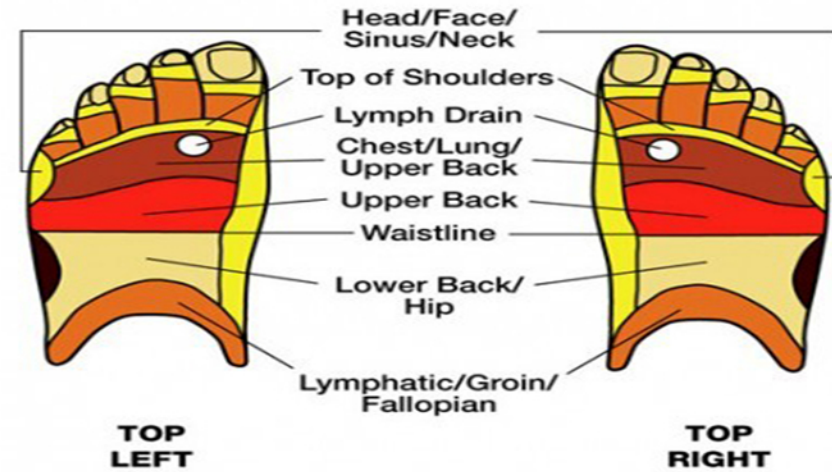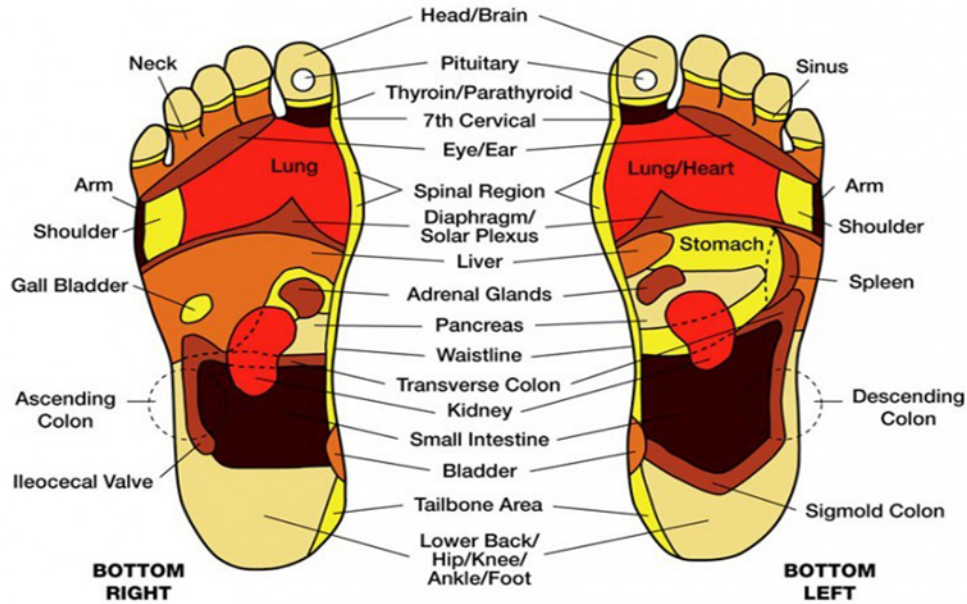


Warming Meridians

The energy of microcrystalline magnetic resonance is similar to that of human body
So as to help the human body warm up from the bottom of the foot to the six meridians
Dissolve stasis, dredge Qi and Blood,
improve immunity

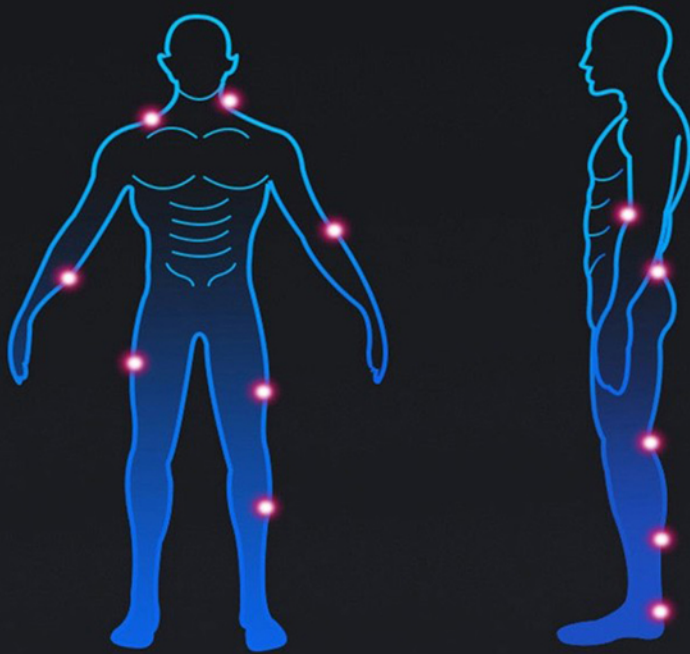

BODY MERIDIANS

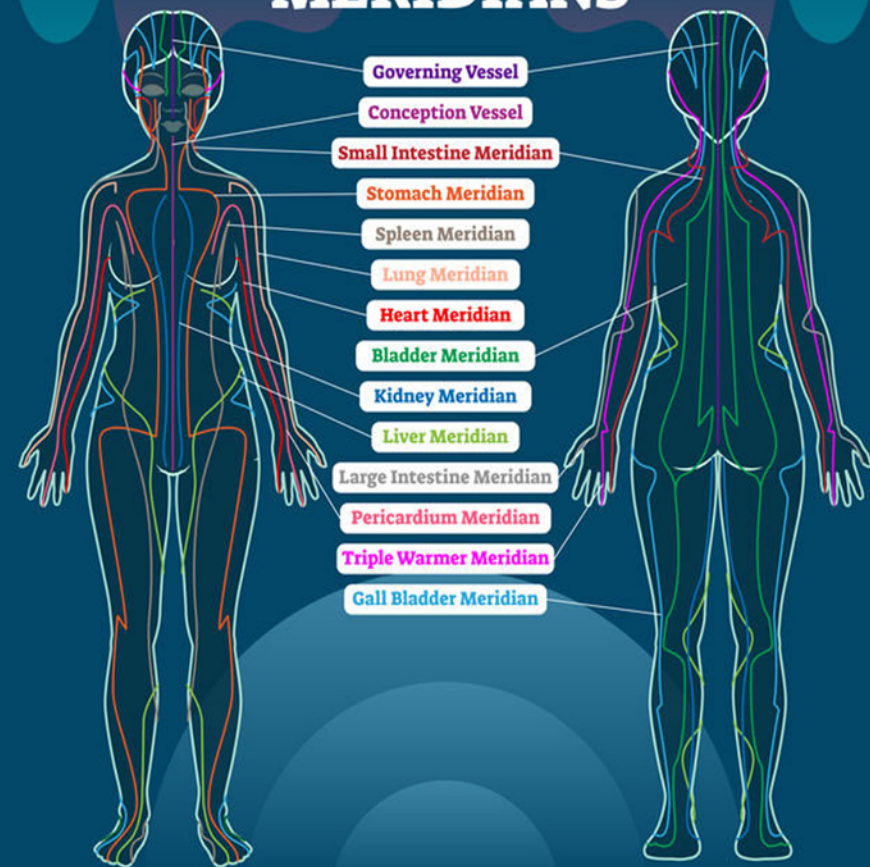


What happens with Terahertz Technology


HOW TO QUICKLY
UNDERSTAND "HOT TECHNOLOGY"



DIRECT INTRODUCTION
OF LIGHT WAVES INTO THE HUMAN BODY
AUTOMATIC IDENTIFICATION

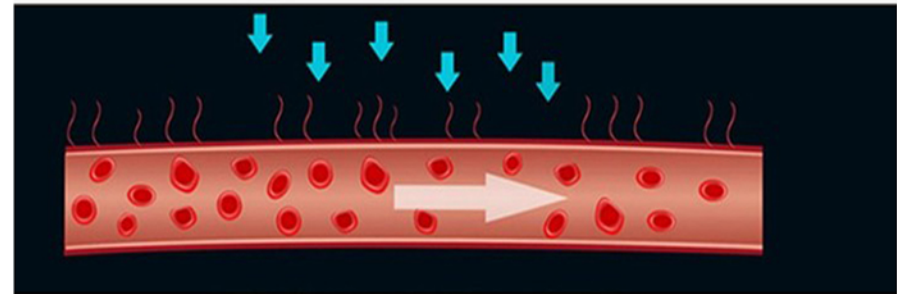

CELLS NORMALLY INCREASE
BLOOD FLOW (NORMAL COOLING)

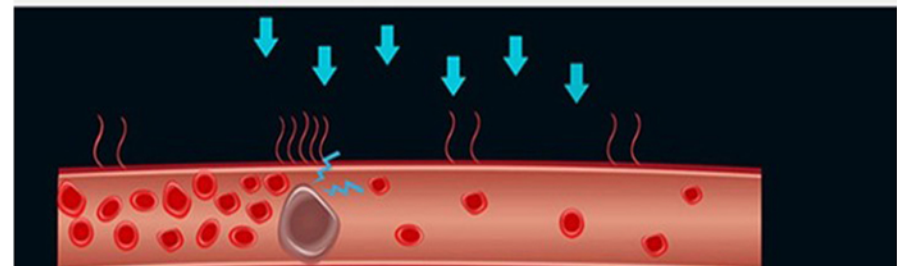

A large, irregular splash of teal and light blue watercolor paint serves as the background for the text. The colors are layered, with darker teal in the center and lighter, almost white, tones at the edges, creating a soft, painterly effect.

Thank You